


MOLDED SILICONE ASSEMBLIES

PLATINUM-CURED SILICONE TUBING AND FITTINGS

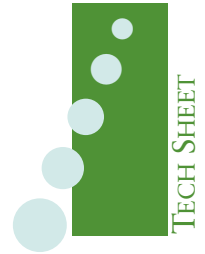

Applied as silicone fittings attached directly to silicone tubing or hose, AdvantaPure's Molded Silicone Assemblies are well-suited for applications involving aseptic transfer, sterile processing, and other critical applications in the pharmaceutical, biopharm, and high purity industries. The assemblies reduce installation time, improve cleanliness, and retain the many benefits of silicone standards. Full BPOG standardized extractables test protocol (six extraction solvents, tested to four different exposure times, out to 70 days) is in process with a fully-accredited, third-party test lab.

Molded Silicone Assemblies are not intended for implantation or continuous steam applications.



KEY FEATURES



- Particulate testing per USP 788 is performed quarterly
- The compound used to make Molded Silicone Assemblies has a Master File with the U.S. Food and Drug Administration
- The assemblies eliminate entrapment issues and leaks associated with barb-type fittings
- Compatible with AdvantaSil™ APST and APHP (high pressure) tubing and with APSH braid-reinforced hose; APHP applications include TFF, virus filtration, inline integrity testing of filtration assemblies, and high volume fluid transfer
- Provides a smooth inner surface with a seamless transition
- Offers consistent interior dimensions for maximum flow rates
- Reduces end-user assembly times—arrives ready to install
- No product contact with metallic materials
- Certified free of animal-derived ingredients
- Retains silicone's general temperature range: -100°F to 400°F
- Extremely flexible
- Sterilizable by CIP, SIP, and autoclave processes
- Imparts no taste or odor to critical streams
- Clean room manufactured
- Sealed in clear polybags during manufacturing process to ensure cleanliness
- Documented lot traceable with identification on bags
- Documented quality control

SILICONE

www.advantapure.com


MOLDED SILICONE ASSEMBLIES

PLATINUM-CURED SILICONE TUBING AND FITTINGS

Molded Silicone Assemblies are available from AdvantaPure in mini and standard Tri-Clamp® fitting styles and are supplied with male gaskets molded directly to the face of the fitting. Straight-through designs are produced in sizes up to 1" tube I.D. Protective backup cups, offered in polysulfone or series 304 stainless steel for extra durability,

provide a stable clamping surface and safeguard the fittings during installation and use.

All assemblies are custom made to best meet each application's requirements and incorporate either unreinforced silicone tubing, including our high pressure style, or braid-reinforced silicone hose. Working pressure ratings vary depending on product. Call AdvantaPure today to discuss your needs.

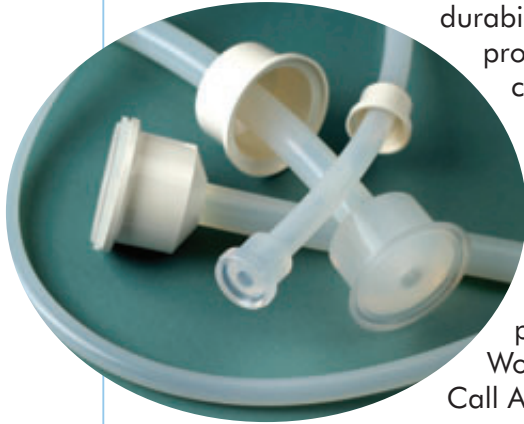

SILICONE

Purity in Fluid Flow Systems®

Molded Silicone Assemblies are part of a complete line of products for the Pharmaceutical, Biomedical, Cosmetic, and Food industries. Items available include . . .



- Unreinforced, 50 Shore A, platinum-cured silicone tubing
- Braid-reinforced, platinum-cured silicone hose
- Wire and mesh reinforced, mandrel-wrapped, platinum-cured silicone suction hose
- AdvantaFlex® biopharmaceutical grade TPE tubing
- Stainless steel overbraided PTFE hose and wire-reinforced EPDM hose
- Sanitary fittings and complete hose assemblies including molded silicone ends
- Silicone injection molding capabilities – manifolds, stoppers, and sealing systems for Single-Use processes

Address: NewAge Industries AdvantaPure®
145 James Way
Southampton, PA 18966 U.S.A.

Toll Free Phone: **888-755-4370**
Phone: **215-526-2151**
Toll Free Fax: **888-258-4293**
Fax: **215-526-2167**
Web Site: **www.advantapure.com**
E-mail: **sales@advantapure.com**

AdvantaPure®
NEWAGE INDUSTRIES



DISTRIBUTED BY:

©2018 NewAge® Industries, Inc. Specifications are subject to change without notice. AdvantaPure products are produced by NewAge® Industries, Inc. AdvantaFlex®, AdvantaSil™, Purity in Fluid Flow Systems®, NewAge Industries AdvantaPure®, NewAge® and the "N" logo are trademarks of NewAge® Industries, Inc. Tri-Clamp® registered trademark of Alfa Laval Inc.

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term NewAge.



ESOP Employee Owned
For Your Benefit