


# BIOREACTOR TUBING KITS

SINGLE-USE & CUSTOMIZABLE FOR BENCHTOP BIOREACTORS

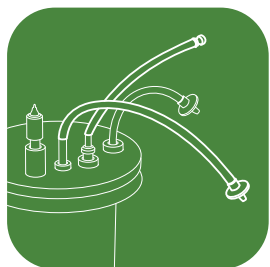

SUS  
www.advantapure.com

Bioreactor Tubing Kits are built-to-order and clean room packaged. The Single-Use kits are custom manufactured and assembled from a variety of components — tubing, filters, connectors — to offer faster turnaround time and increased processing capacity. Made from platinum-cured silicone or AdvantaFlex® biopharmaceutical grade TPE tubing, AdvantaPure's Bioreactor Tubing Kits provide convenience and greater efficiency.


Full BPOG standardized extractables test protocol has been performed by a fully-accredited, third-party test lab; applies to tubing only.

## KEY FEATURES

- Components meet relevant ISO and USP standards, including Class VI
- Assembled and packaged in ISO 7 clean rooms
- AdvantaFlex TPE tubing is sterile weldable and heat sealable
- Tubing is certified free of animal-derived ingredients
- Documented lot traceable with identification on packaging
- Documented quality control per ISO 9001:2015
- Tubing material validation and extractables test portfolios available upon request
- Optional RFID tags easily attach to assemblies for identification purposes

# BIOREACTOR TUBING KITS

SINGLE-USE & CUSTOMIZABLE FOR BENCHTOP BIOREACTORS

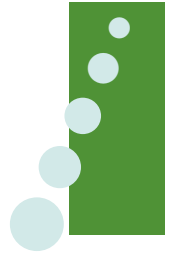

## AdvantaPure's Bioreactor Tubing Kits:



- Offer faster turnaround time for your bioreactor and increase your processing capacity
- Decrease cross contamination risks
- Reduce repetitive motion injuries from assembly processes
- Minimize the number of inventoried components you need to stock
- Provide a single, animal-derived-ingredient-free material contact surface (both tubing and molded components)
- Design your ready-to-use assembly to your specifications


SUS

■ Kits with AdvantaFlex® TPE tubing allow for easy welding and sealing.



## Bioreactor Tubing Kits are Custom Assembled for Your Process

- Specify your tubing material and length for addition, vent, sparge and harvest lines.

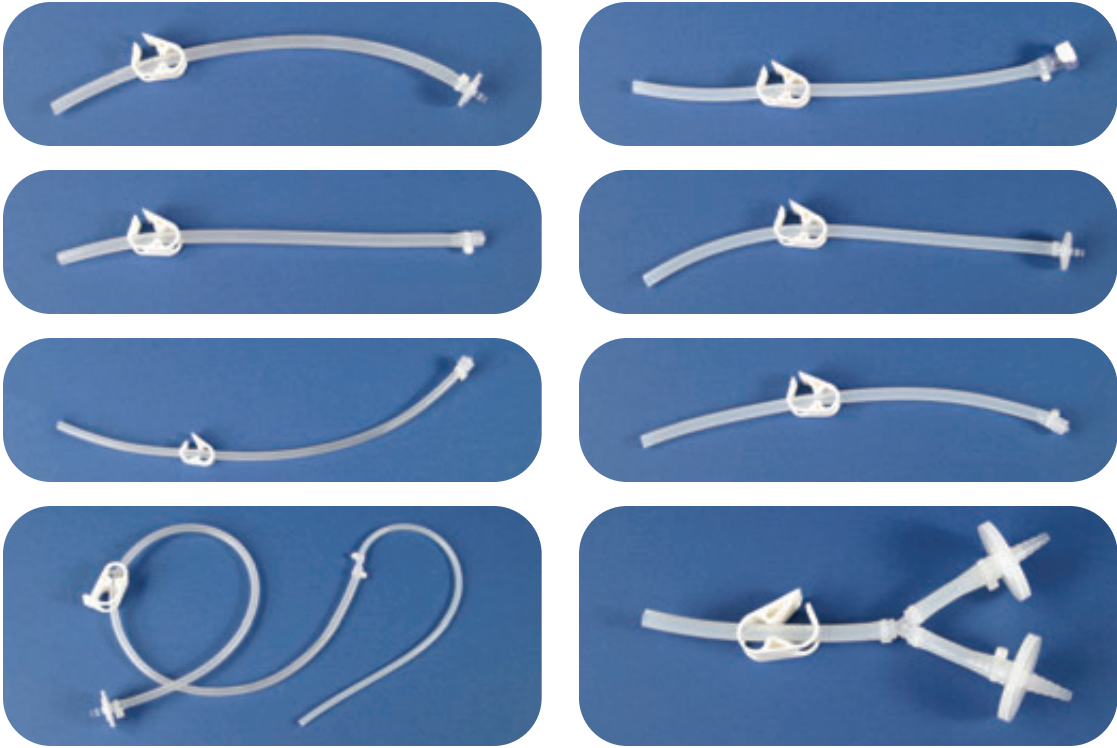

- You identify the tubing sizes, materials and lengths, plus filters and connectors.

SUS




- AdvantaFlex® TPE lines are available with sealed ends; silicone lines are available plugged.




# BIOREACTOR TUBING KITS

SINGLE-USE & CUSTOMIZABLE FOR BENCHTOP BIOREACTORS

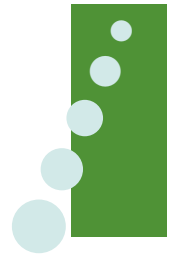

Select from AdvantaFlex TPE Tubing or APST Silicone Tubing

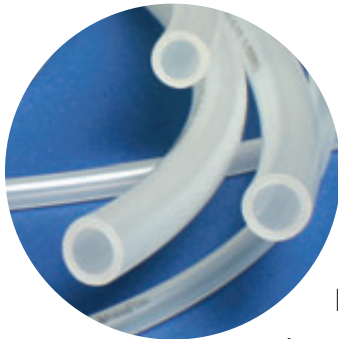

 **AdvantaFlex**<sup>®</sup>  
BioPharmaceutical Grade Tubing



AdvantaFlex addresses the need for a flexible, sterilizable, moldable, weldable and heat sealable biopharmaceutical tubing without silicone oils. AdvantaFlex maintains its physical properties following sterilization processes, remains translucent for visible product flow and does not become gummy. The tubing resists kinking and is an excellent component for Single-Use systems.

Learn more at [www.advantapure.com/advantaflex.htm](http://www.advantapure.com/advantaflex.htm)



## APST UNREINFORCED SILICONE TUBING

Low-volatile grade, platinum-cured silicone tubing is cleanroom produced for critical pharmaceutical applications. APST has undergone extensive physical, chemical, and biological testing and meets USP Class VI, FDA CFR 177.2600, ISO 10993, and European Pharmacopoeia 3.1.9 standards.



Learn more at [www.advantapure.com/apst-silicone-tubing.htm](http://www.advantapure.com/apst-silicone-tubing.htm)

SUS

Address: NewAge Industries AdvantaPure<sup>®</sup>  
145 James Way  
Southampton, PA 18966 U.S.A.

 **AdvantaPure**<sup>®</sup>  
NEWAGE INDUSTRIES



Toll Free Phone: **888-755-4370**  
Phone: **215-526-2151**  
Toll Free Fax: **888-258-4293**  
Fax: **215-526-2167**  
Web Site: **[www.advantapure.com](http://www.advantapure.com)**  
E-mail: **[sales@advantapure.com](mailto:sales@advantapure.com)**

DISTRIBUTED BY:



**ESOP** Employee Owned  
for Your Benefit

©2021 NewAge<sup>®</sup> Industries, Inc. Specifications are subject to change without notice. AdvantaPure products are produced by NewAge<sup>®</sup> Industries, Inc. AdvantaFlex<sup>®</sup>, Purity in Fluid Flow Systems<sup>®</sup>, NewAge<sup>®</sup>, NewAge Industries AdvantaPure<sup>®</sup>, and the "N" logo are registered trademarks of NewAge<sup>®</sup> Industries, Inc.

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term New-Age.